


Hot Backup is a cost-effective solution that enables users to actively mirror the Instrument Manager™ database. It provides real time, continuous data backup, 24/7, without user intervention.



In the Event of PC Failure... users can follow the comprehensive step-by-step instructions to bring the shadow system into live operation in less than ten minutes. The new PC will take over exactly where the primary system left off.

How It Works... On the Primary system, transactions are stored in journal files and then applied to the Primary system's Caché database i.e., its cache.dat files.

On the Shadow system journal files are copied from the Primary system where those copied journal files' are applied to the Shadow system's Caché database.

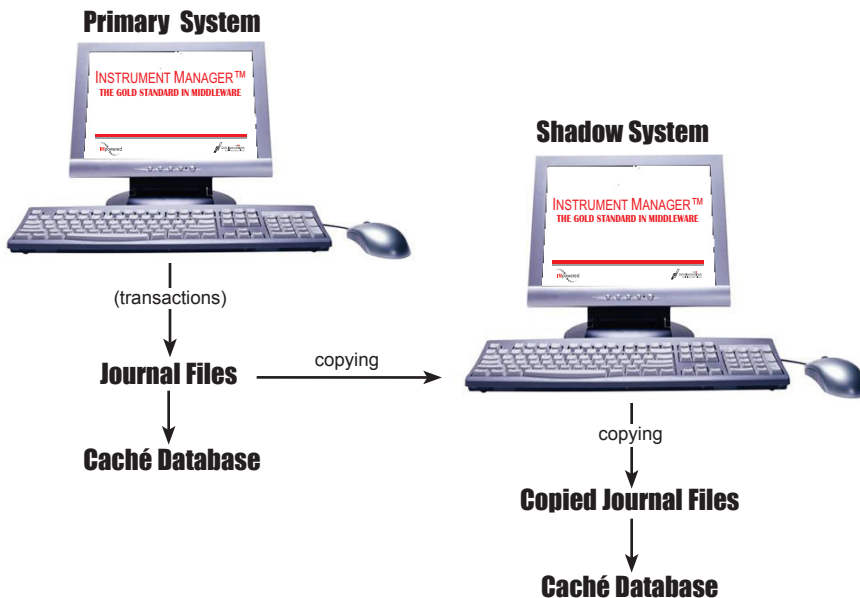

Fast Facts: Hot Backup...

- Actively mirrors data
- Minimal downtime to failover
- Cost effective
- Comprehensive documentation

Data Innovations offers the most comprehensive Middleware system in the market to manage laboratory operations, including pre-analytical, analytical, and post-analytical sample processing and non clinical tasks such as equipment maintenance and specimen archiving. Instrument Manager™ (IM), the Gold Standard in Middleware, is available independent of instrument or LIS vendor, laboratory department, and clinical discipline.

DATA INNOVATIONS, INC.®

Phone: (802) 658-1955
Fax: (802) 658-2782
Email: sales@datainnovations.com

DATA INNOVATIONS EUROPE

Phone: 32 (0)2 332 24 13
Fax: 32 (0)2 376 43 84
Email: europe-sales@datainnovations.com

DATA INNOVATIONS ASIA

Phone: (852) 2398 3182
Fax: (852) 2398 8667
Email: asia-sales@datainnovations.com

DATA INNOVATIONS LATIN AMERICA

Phone: 55-11-3801-3283
Fax: 55-11-3871-9592
Email: latinamerica-sales@datainnovations.com