

# Integrated QC Solution


Providing real-time integration between  
QC and results production

The reputation of your lab depends on the integrity of your results. At Data Innovations we recognize the importance of integrating QC into results production. Nowhere is this more critical to your business than in quality control procedures that protect your physicians and their patients.

Data Innovations is the only middleware company to offer a bi-directional interface with Bio-Rad's QC OnCall™ and Unity Real Time®.

## Other features of Data Innovations' Integrated QC include:

- Automatically and immediately sends QC results to QC OnCall™ / Unity Real Time®
- Alter results production immediately based on Bio Rad's Unity Real Time® evaluation and user-defined QC rules in Instrument Manager™
- On-demand interlaboratory comparison reports for real-time troubleshooting
- Web-based, peer grouping via Bio-Rad's QCNet™
- Input QC comments into Unity Real Time® directly from Instrument Manager users screens

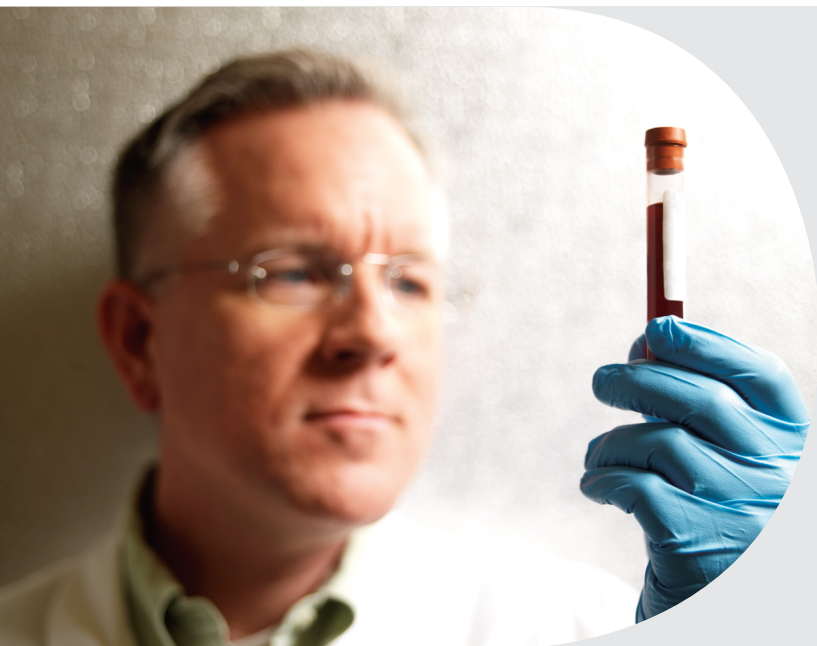

## Benefits of Notifier

Along with the obvious benefits of Integrated QC is Data Innovations' unique, powerful Notifier application. Be alerted immediately with pop-up notifications that deliver critically flagged information to operators, for example, when a 1-3s rule violation occurs or QC has not been run in 8 hours for a particular analyte, so they can take appropriate action to maintain your quality standards.